

BIOMÉRIEUX

VIDAS® KUBE™

Enter a World of Simplicity

TRUSTED. ROBUST. UNIVERSAL.



PIONEERING DIAGNOSTICS

VIDAS® KUBE™

YOUR IMMUNOASSAY SOLUTION

Your lab immunoassay needs are evolving. You need to provide same-day, accurate results in a timely manner to ensure patient and clinician satisfaction remains high.

Enter the world of simplicity with VIDAS® KUBE™ a truly flexible, cost-effective automated platform that provides rapid results in complete confidence. Scaling with ease, VIDAS® KUBE™ is a simple and reliable addition to your evolving lab!



EASY IMPLEMENTATION

Adding VIDAS® KUBE™ to your lab is simple, benefit from:

- Accessible volume flexibility - up to 36 tests/hour per module
- Stackable benchtop design - 1-6 modules
- Minimum training & installation required
- Integrated computer for up to 2 modules

INCREASED CONFIDENCE IN RESULTS

VIDAS® KUBE™ ensures your lab has reliable and accurate results every time.

- On-board Strip & SPR compatibility verification
- Optical sample detection
- Alerts in case of sample issues (leakage/clogging)
- Shared calibration between modules

VIDAS® KUBE™ facilitates easy user adoption through its:

- Intuitive & easy to use software interface
- Bi-directional LIS² connectivity
- Secure data with cybersecurity protection
- Full traceability (21CFR³, Audit trail)

RECOGNIZED HIGH QUALITY SERVICE

VIDAS® KUBE™ offers long-term peace of mind with:

- Evolutive and adaptable training journey
- Optimized productivity with remote maintenance via VILINK® connectivity
- Reduced maintenance and controls (no need for QCV®¹)


VIDAS® KUBE™ adaptable automation will help your lab be more productive, more accurate, and more innovative - all in one simplified platform.



VIDAS® CONCEPT INSIDE


UNIQUE SINGLE TEST CONCEPT

- Ready to use reagents with long shelf-life
- Everything is included in the kit

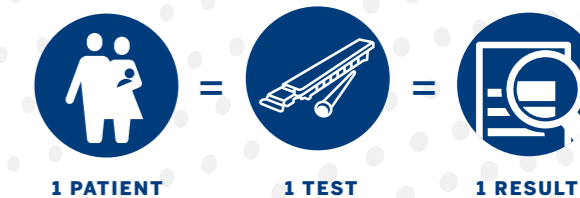

BROAD ROUTINE TEST MENU WITH MORE THAN 80 ASSAYS

- Performance demonstrated in analytical studies*
- VIDAS® Assay Solutions cover:
 - Emergency & Critical Care
 - Immunochemistry
 - Infectious Diseases

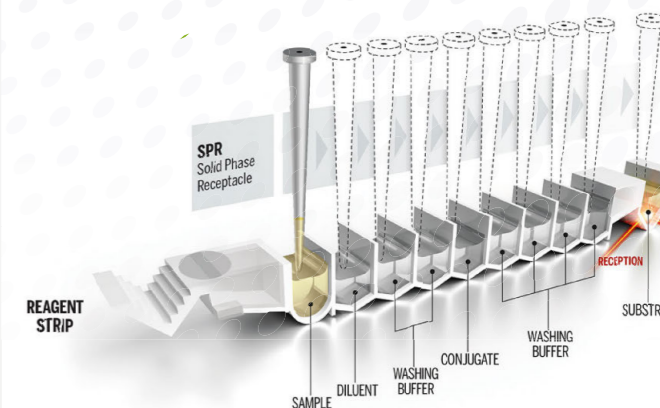

¹ Quality Control VIDAS®
² Laboratory Information System
³ Code of Federal Regulations

*Please consult your local representative to find out more

THE WORLD OF VIDAS®


We understand your **evolving lab needs**, whether you are new to automation, or upgrading to advanced technology.

We recognize your need for **simple, built-to-last** solutions with excellent service to keep your lab up and running.

We designed VIDAS® KUBE™ so you could **benefit from advanced technology** while keeping everything labs appreciate about the recognized VIDAS® solution.

By partnering with bioMérieux, you gain **over 30 years** of worldwide experience in immunoassays.

**OVER 3 DECADES
OF INNOVATION**

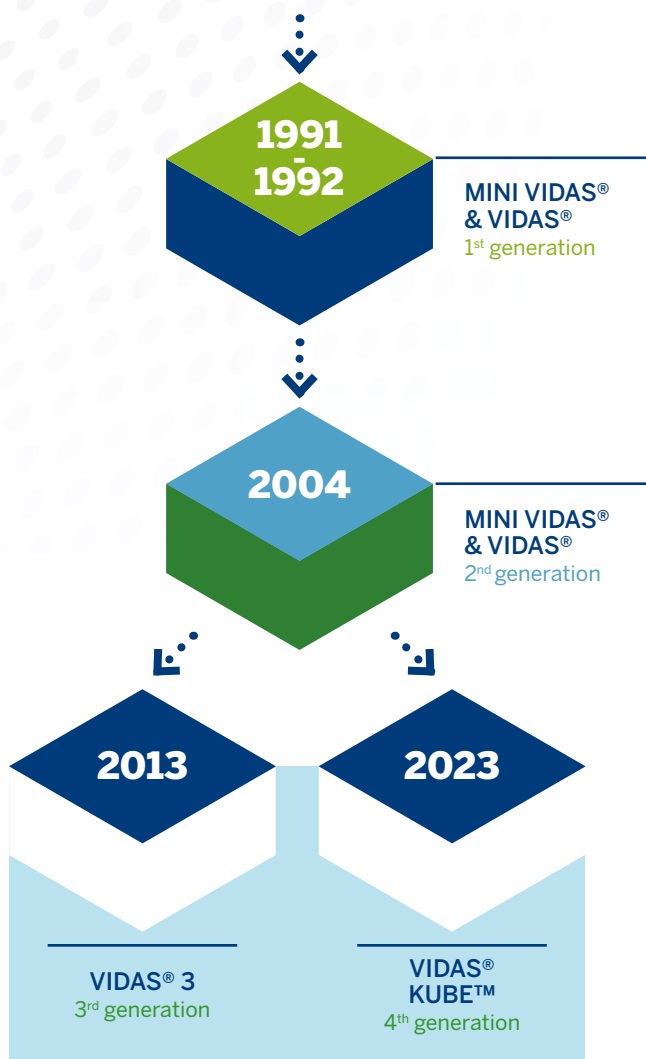

INSTRUMENT SPECIFICATIONS

VIDAS® KUBE™ References	423 912 (single unit) 470 7576 (standard configuration)
Weight	23.7 kg (single unit)
Dimensions	L. 43.4 cm x H 32.1 cm x P 48.5 cm (single unit)
Electrical	100/240 V, 56-60Hz
Power	150 W
Energy efficient	Sleep mode; Programmable morning start up
Operating temperature	18-30°C
Operating humidity	20-80% relative humidity non-condensing
Altitude maximum	2700 m
Sound level in operating mode	≤ 70 dBa
Usability criteria	<ul style="list-style-type: none"> • Data purge: Automatic data archiving • User safety during testing: Auto-locking of section doors • Calibration setting: Modifiable

PC specifications

An additional PC is required for the VIDAS® KUBE™ configuration from 3 to 6 modules.

OS Configuration	Windows 10
CPU	Intel Core i3 – 8100 3.6 GHz
Memory	27000 tests results Approx. 5 years of data
RAM	DDR4 - 16GB 2400 MHz
Connectors	USB port x 6 Ethernet port x 1 Display port (screen) x 1

10/2022 / 9322986.008/GB/A / This document is not legally binding. bioMérieux reserves the right to modify specifications without notice. / BIOMÉRIEUX, the BIOMÉRIEUX logo, PIONEERING DIAGNOSTICS, VIDAS, VIDAS KUBE, VILINK are used, pending and/or registered trademarks belonging to bioMérieux or one of its subsidiaries or one of its companies. Any other name or trademark is the property of its respective owner / Photos : Getty Images / bioMérieux S.A. RCS Lyon 673 620 399 / Printed in France /thera/RCS Lyon B 398 160 242 : Do not litter